


# Application note #7090


## Separation of Psoralen and Isopsoralen

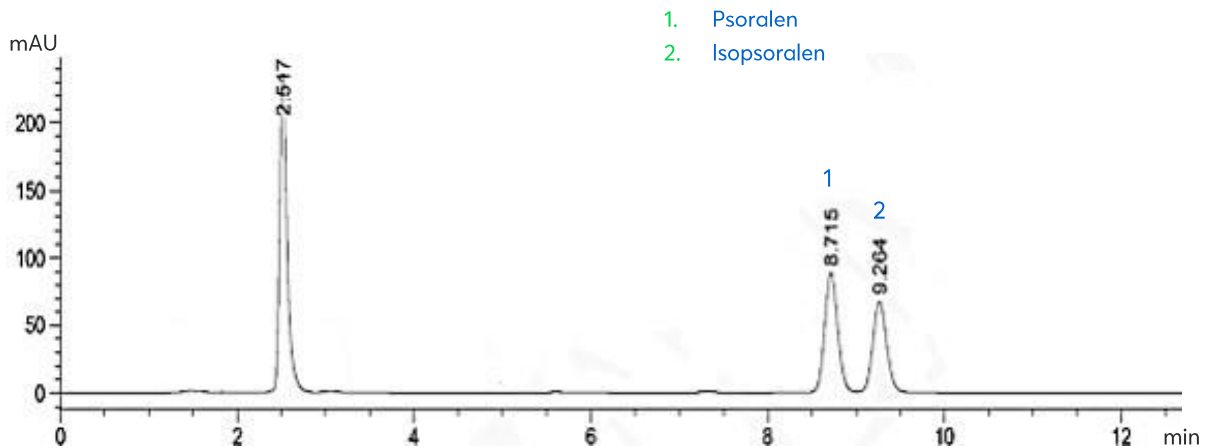


# Method Details

## CONDITIONS

Column:	Avantor® ACE® Excel 5 C18-PFP
Particle Size:	5 µm
Dimensions:	250 x 4.6 mm
Mobile Phase:	H <sub>2</sub> O/MeOH (40:60 v/v)
Flow Rate:	1.0 mL/min
Injection:	10 µL
Detection:	UV, 240 nm

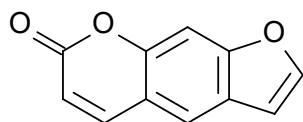

Psoralen

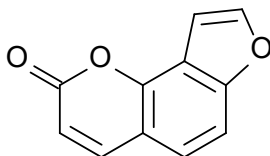

Isopsoralen

## ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® Excel 5 C18-PFP	HPLC Column	250 x 4.6 mm	EXL-1210-2546U