
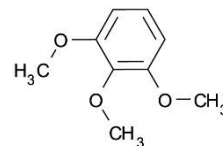



Application #AN2430

Conditions


Column: ACE 3 C18-PFP
Dimensions: 150 x 4.6 mm
Part Number: ACE-1110-1546
Mobile Phase: H₂O/MeOH (50:50 v/v)
Flow Rate: 1 mL/min
Injection: 5 µL
Temperature: 40 °C
Detection: UV, 214 nm

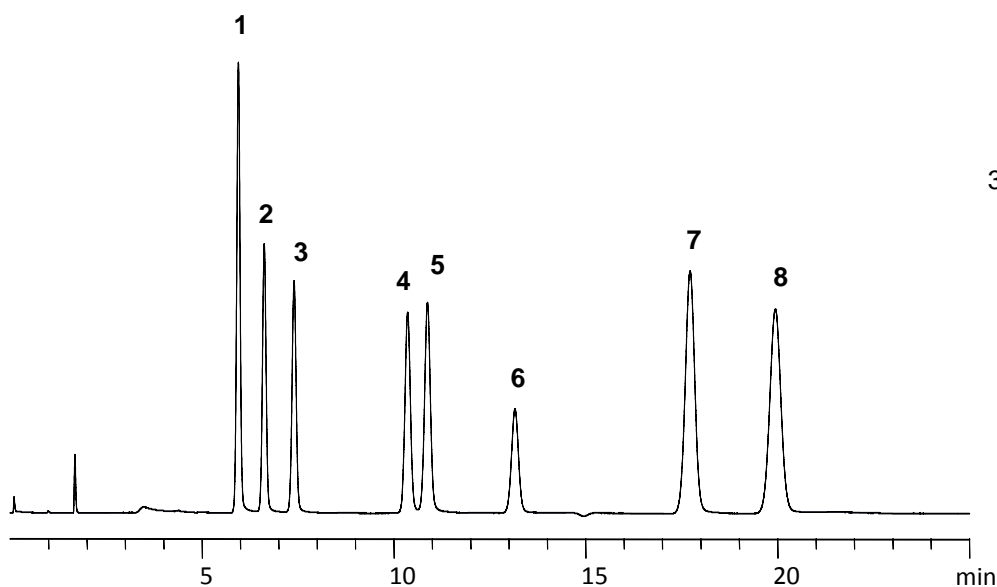


1. 1,2,3-Trimethoxybenzene

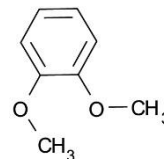


2. 1,2-Dimethoxybenzene

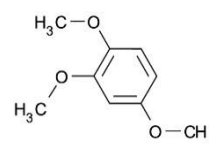


3. 1,2,4-Trimethoxybenzene

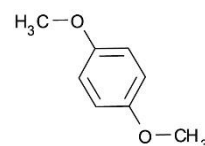


4. 1,4-Dimethoxybenzene

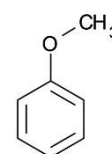


5. Methoxybenzene

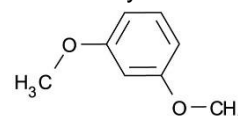

6. 1,3-Dimethoxybenzene

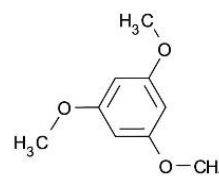

7. 1,3,5-Trimethoxybenzene

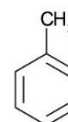

8. Toluene