



Conditions


Column:	ACE Excel 5 NH2
Dimensions:	250 x 4.6 mm
Part Number:	EXL-1214-2546U
Mobile Phase:	MeCN:H ₂ O 80:20 v/v
Flow Rate:	1.5 mL/min
Temperature:	35 °C
Detection:	Refractive Index
Sample:	Carbo Energy Drink and Sugar standards

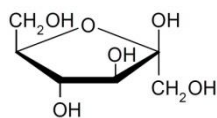


1. Fructose

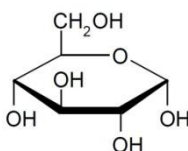


2. Glucose

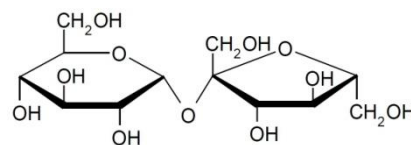


3. Sucrose

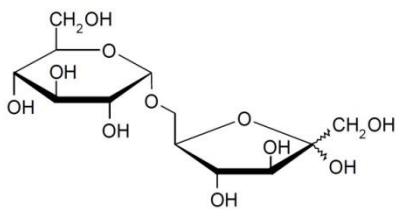

4. Palatinose

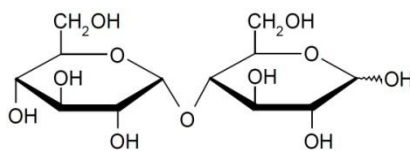


5. Maltose

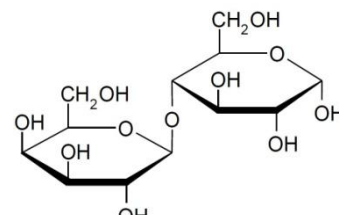


6. Lactose

Blue Exotic Carbo Energy Drink

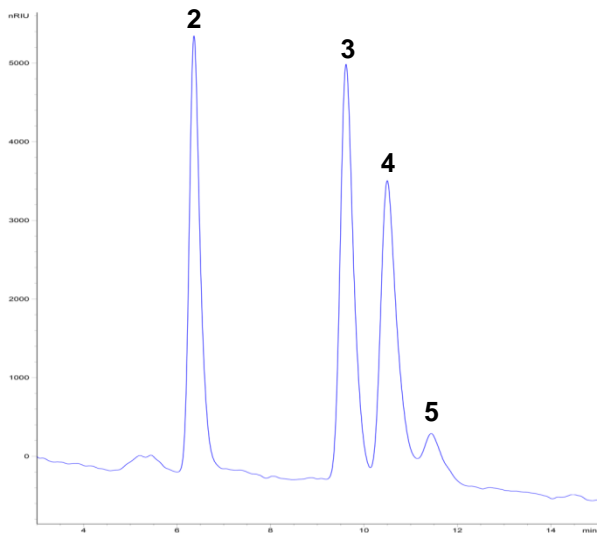

Sugar standards

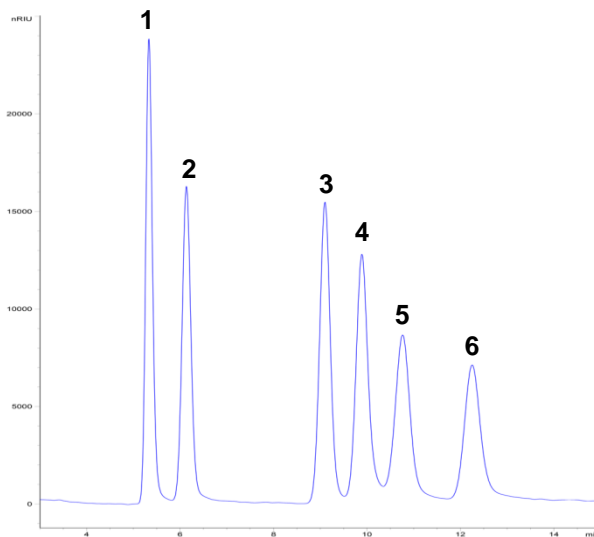

Reproduced with permission of Clinical Nutrition, S.A.U., Barcelona, Spain

© 2016 ACT Ltd. All rights reserved. All products are available worldwide.