


Conditions

Column: ACE Excel 3 C18-PFP
 Dimensions: 150 x 4.6 mm
 Part Number: EXL-1110-1546U
 Mobile Phase: A: H₂O
 B: MeCN
 C: 10% formic acid

Time (mins)	%A	%B	%C
0	84	15	1
12	74	25	1
14	59	40	1
16	84	15	1
18	84	15	1

Flow Rate: 1 mL/min
 Detection: UV, 268 nm

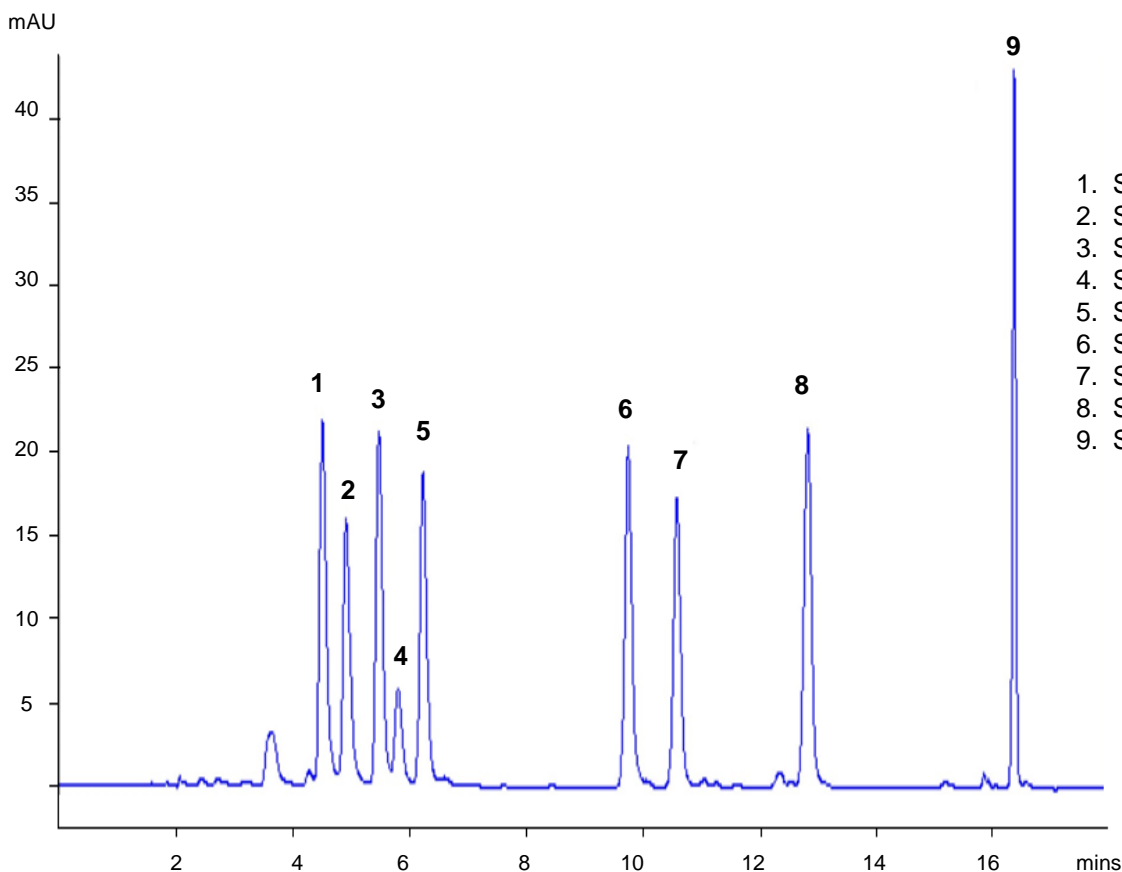

1. Sulfadiazine
2. Sulfapyridine
3. Sulfamerazine
4. Sulfamoxole
5. Sulfamethazine
6. Sulfamonomethoxine
7. Sulfachloropyridazine
8. Sulfamethoxazole
9. Sulfadimethoxine

Reproduced with permission of National Research Institute of Agriculture and Food Technology INIA, Spain

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.