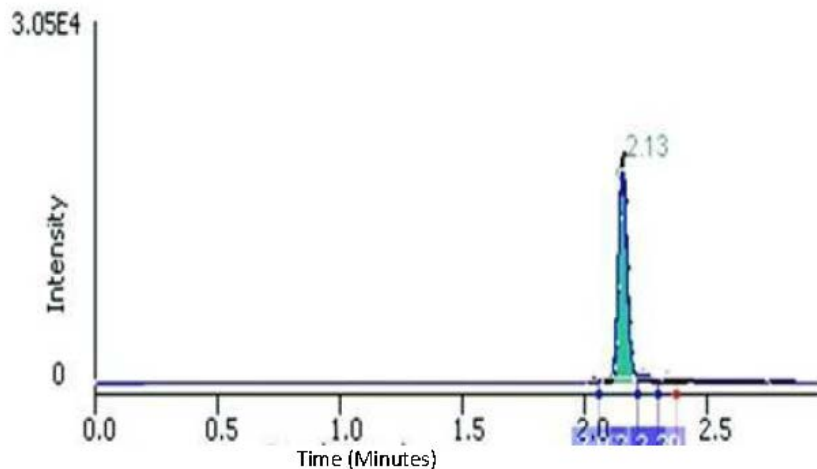
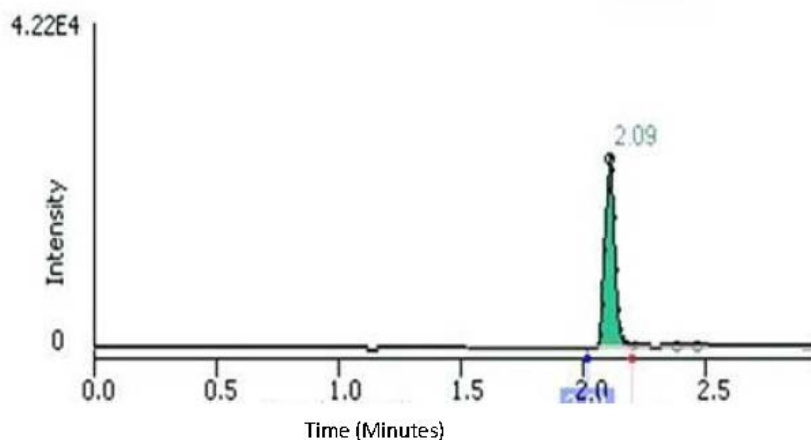




# Testosterones



**Figure 1(A)** LC-MS/MS EIC of testosterone in protein precipitated human serum,



**Figure 1 (B)** LC-MS/MS EIC of 17 -hydroxyprogesterone in protein precipitated human serum.

**Fortis C18 50x2.1mm 3um**

**0.1% formic acid + ammonium acetate in Water**

**0.1% formic acid + ammonium acetate in MeOH**

**0.45ml/min**

**40oC**

**Testosterone**

**17-Hydroxyprogesterone**