


Tocopherols

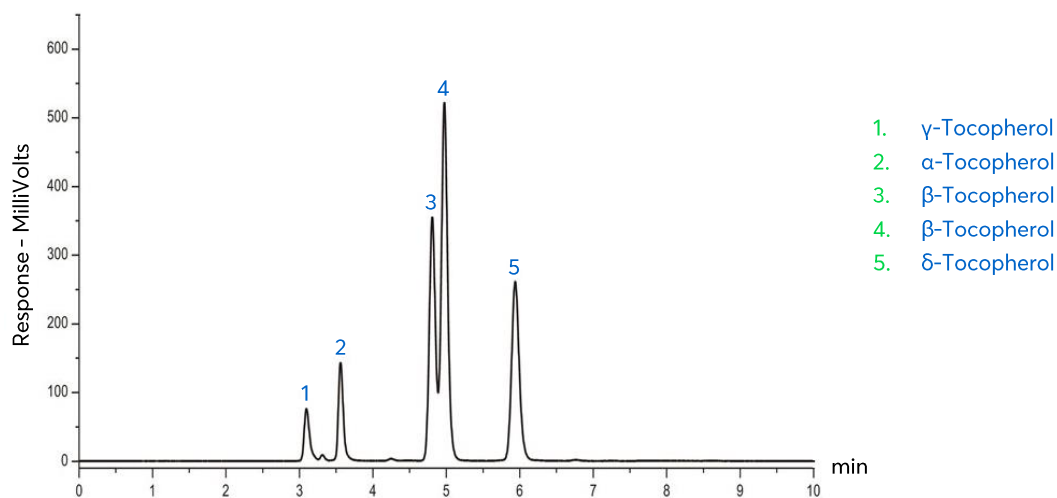

Method Details

CONDITIONS

Column:	Avantor® ACE® 5 SIL
Particle Size:	5 µm
Dimensions:	250 x 4.6 mm
Mobile Phase:	Hexane/IPA (98:2 v/v)
Flow Rate:	1 mL/min
Injection:	1 µL
Temperature:	Ambient
Detection:	UV-Vis, 450 nm

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® 5 SIL	HPLC Column	250 x 4.6 mm	ACE-127-2546
n-Hexane	VWR HiPerSolv CHROMANORM® for HPLC	2.5 L	24575.320P
2-Propanol	VWR HiPerSolv CHROMANORM® for HPLC	2.5 L	20880.320P