
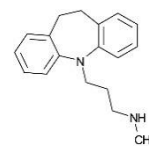


Tricyclic Antidepressants (Isocratic Rapid Analysis)


Application #AN1700

Conditions


Column: ACE 5 C18
Dimensions: 250 x 4.6 mm
Part Number: ACE-121-2546
Mobile Phase: 20 mM ammonium acetate pH 6.3 MeCN/H₂O (65:35 v/v)
Flow Rate: 1.5 mL/min
Injection: 10 µL
Temperature: 60 °C
Detection: UV, 215 nm

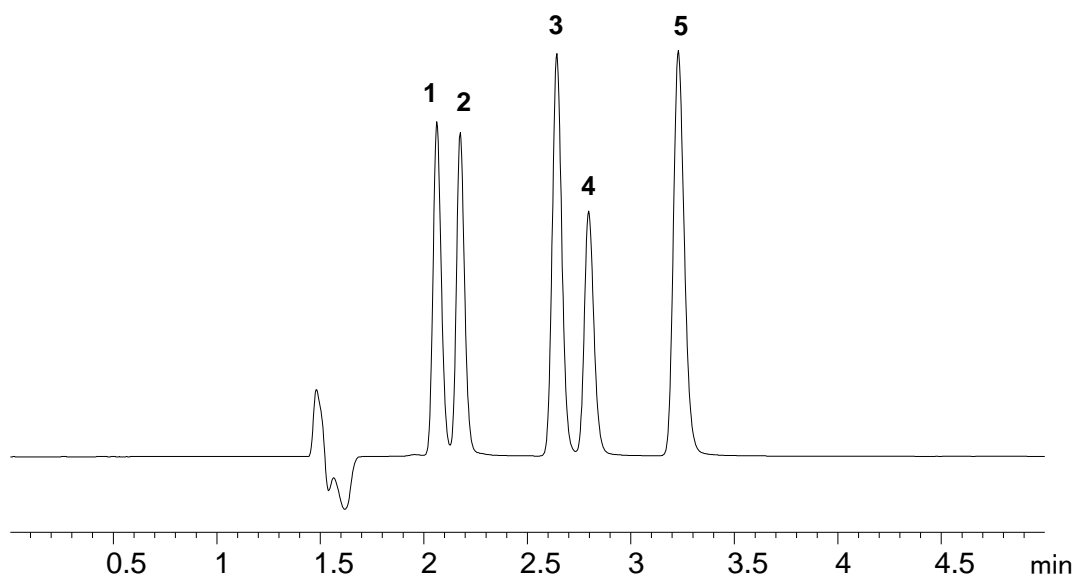


1. Desipramine

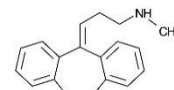


2. Nortriptyline

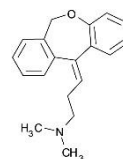

3. Doxepin

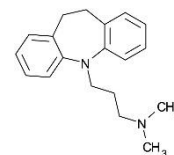

4. Imipramine

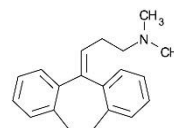

5. Amitriptyline