


# Tricyclic Antidepressants

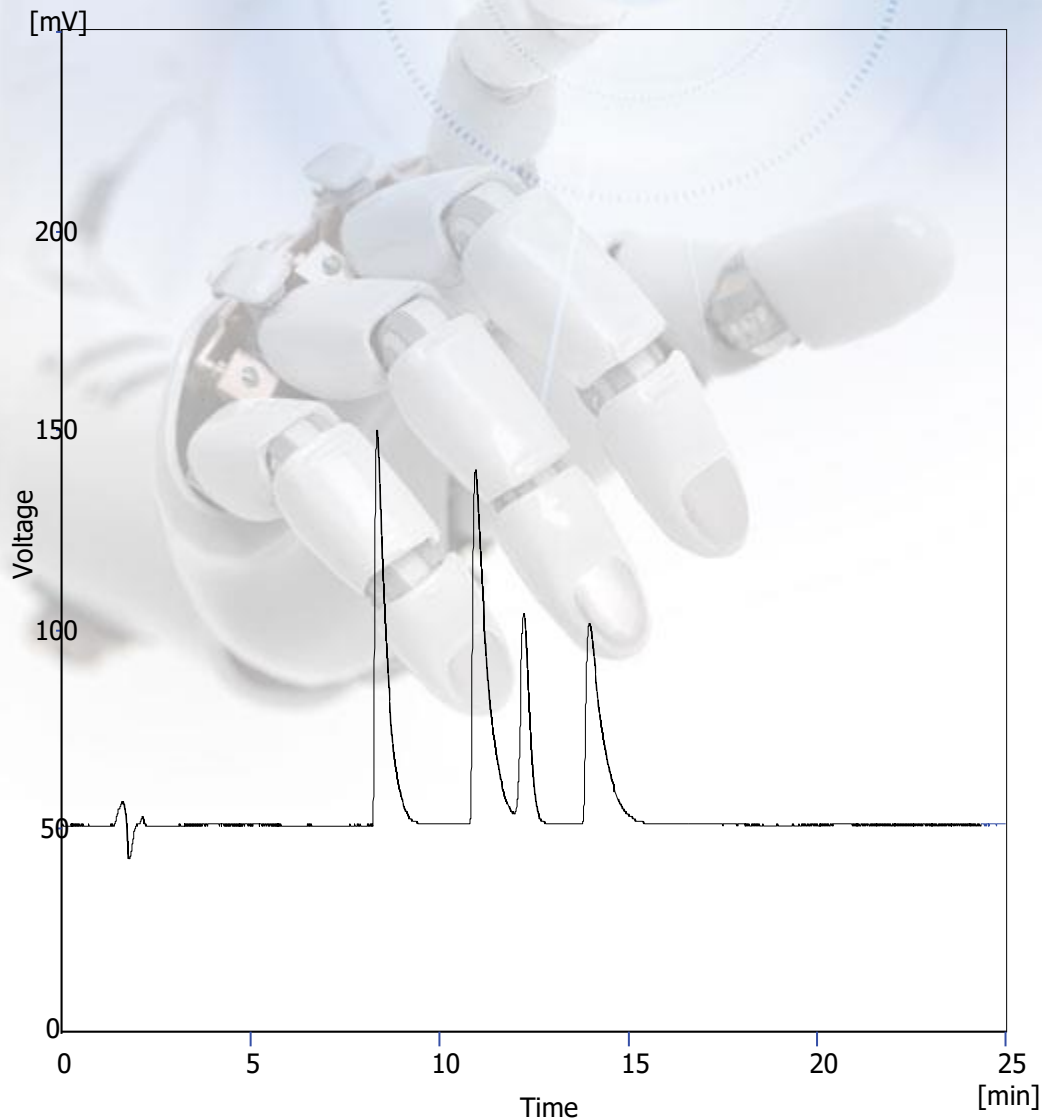

- Column : **Evosphere C18/AR 150x4.6 5um**
- p/n : **EVO18AR-050705**
- Mobile Phase:  
**50:50 KH<sub>2</sub>PO<sub>4</sub> pH 7.0:ACN**

- Flow : **1.0 ml/min**
- Temp : **40°C**
- $\lambda$  : **254nm**

1. **Imipramine Hydrochloride**
2. **Amitriptyline Hydrochloride**
3. **Clomipramine Hydrochloride**
4. **Mianserin Hydrochloride**

