


Tricyclic Antidepressants



- Column : **Evosphere RP18-Amide 150x4.6 3um**
 - p/n : **EVORP18-050703**
 - Mobile Phase:
 - A- 0.1% Formic acid in Water**
 - B- 0.1% Formic acid in ACN**
 - 25-50% B in 10 mins**
 - Flow : **1.0 ml/min**
 - Temp : **30°C**
 - λ : **254nm**
1. **Doxepin**
 2. **Imipramine**
 3. **Nortriptyline**
 4. **Trimipramine**

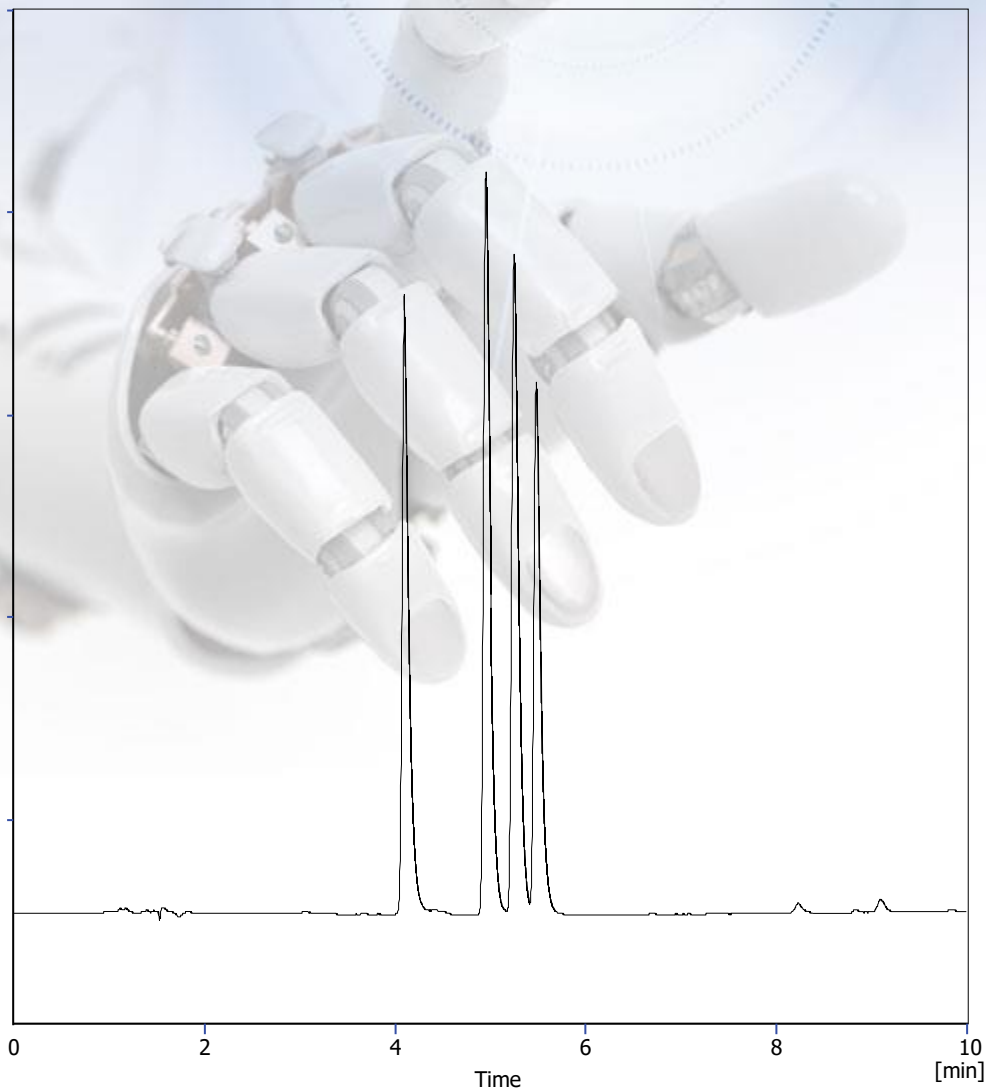