


# USP Method for Analysis of Povidone Impurity A (Vinylpyrrolidinone)

Application #AN6830

## Conditions

Columns:	ACE 5 C18
Dimensions:	150 x 4.6 mm
Part Number:	ACE-121-1546
Mobile Phase:	10:90 MeCN/H <sub>2</sub> O
Flow Rate:	1.0 mL/min
Injection:	5 µL
Temperature:	40 °C
Detection:	UV, 235 nm

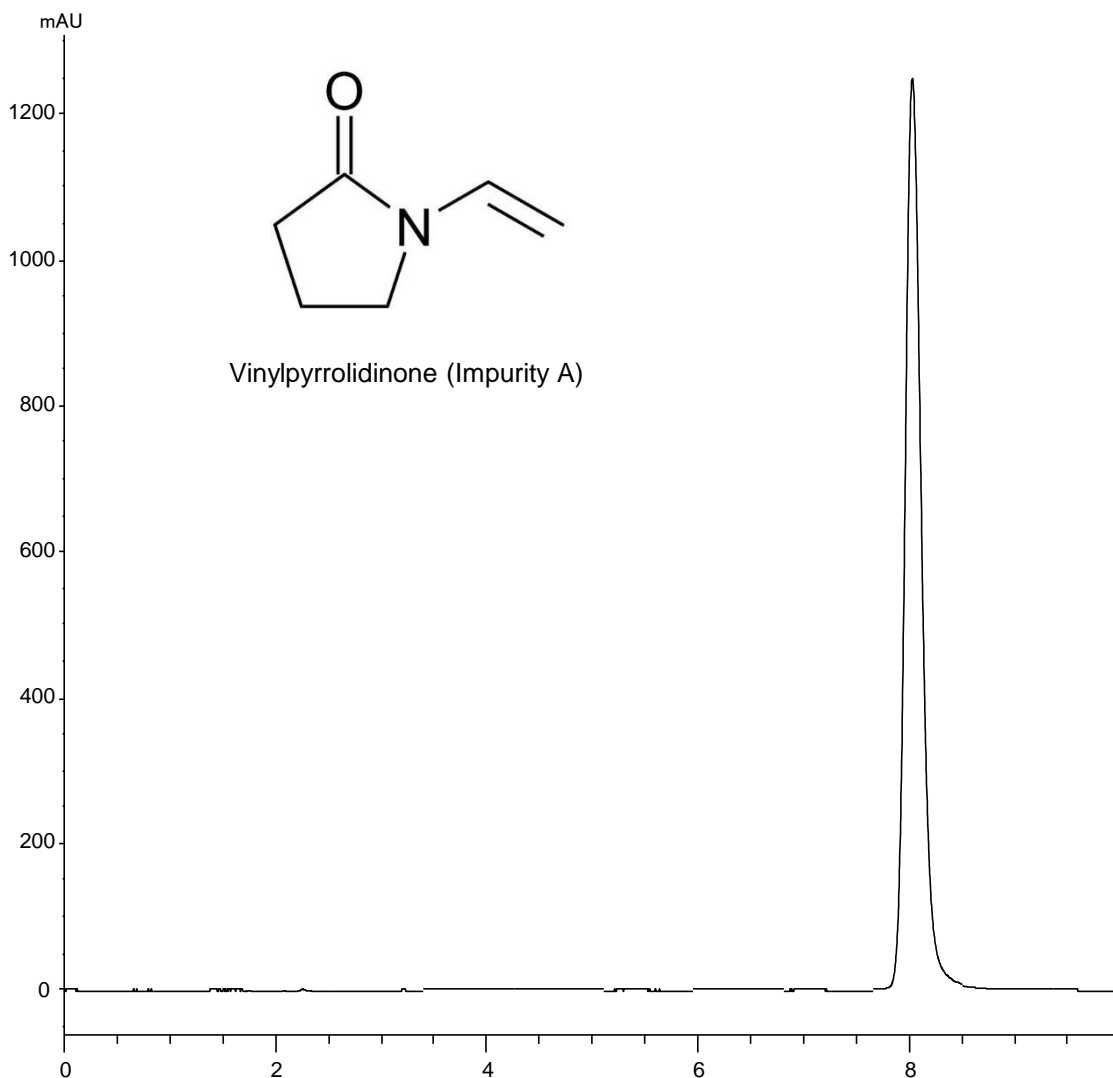

Always check the latest and official method information from the relevant pharmacopoeia prior to analysis.

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit  
[www.ace-hplc.com](http://www.ace-hplc.com) or email: [info@ace-hplc.com](mailto:info@ace-hplc.com)