
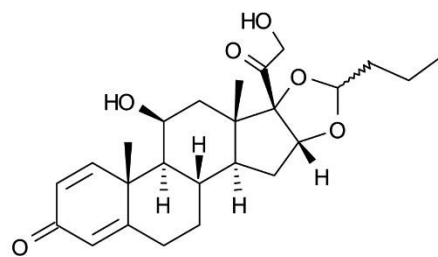


Conditions

Column:	ACE 5 C18
Dimensions:	150 x 4.6 mm
Part Number:	ACE-121-1546
Mobile Phase:	18 mM monobasic sodium phosphate pH 3.2 H ₂ O/MeCN (68:32 v/v)
Flow Rate:	1.5 mL/min
Injection:	20 µL
Temperature:	Ambient (22 °C)
Detection:	UV, 254 nm

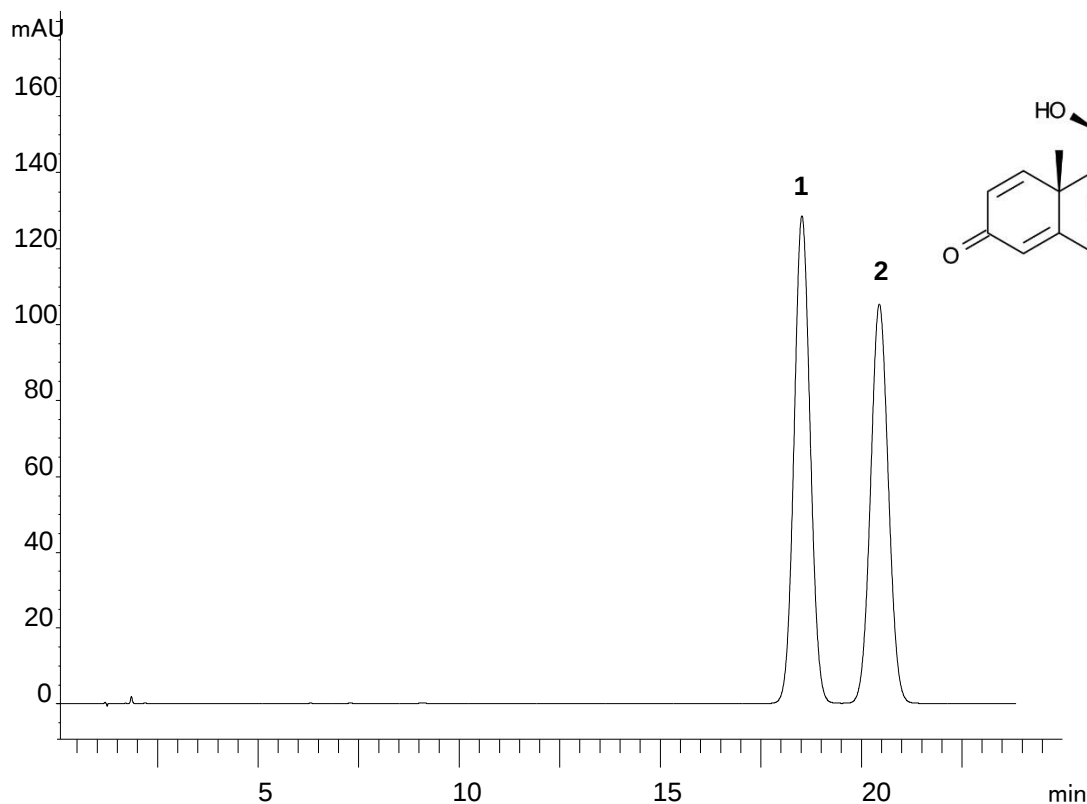

1. Budesonide B
2. Budesonide A

Translated method according to USP <621> guidance. Always check the latest and official method information from the relevant pharmacopoeia prior to analysis.

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.