


Conditions

| | |
|---------------|--|
| Column: | ACE 3 C8 |
| Dimensions: | 50 x 4.6 mm |
| Part Number: | ACE-112-0546 |
| Mobile Phase: | 0.2 M monobasic NaH ₂ PO ₄ pH 2.5 in H ₂ O/MeOH (30:70 v/v) |
| Flow Rate: | 0.56 mL/min |
| Injection: | 4.5 µL |
| Temperature: | 50 °C |
| Detection: | UV, 254 nm |

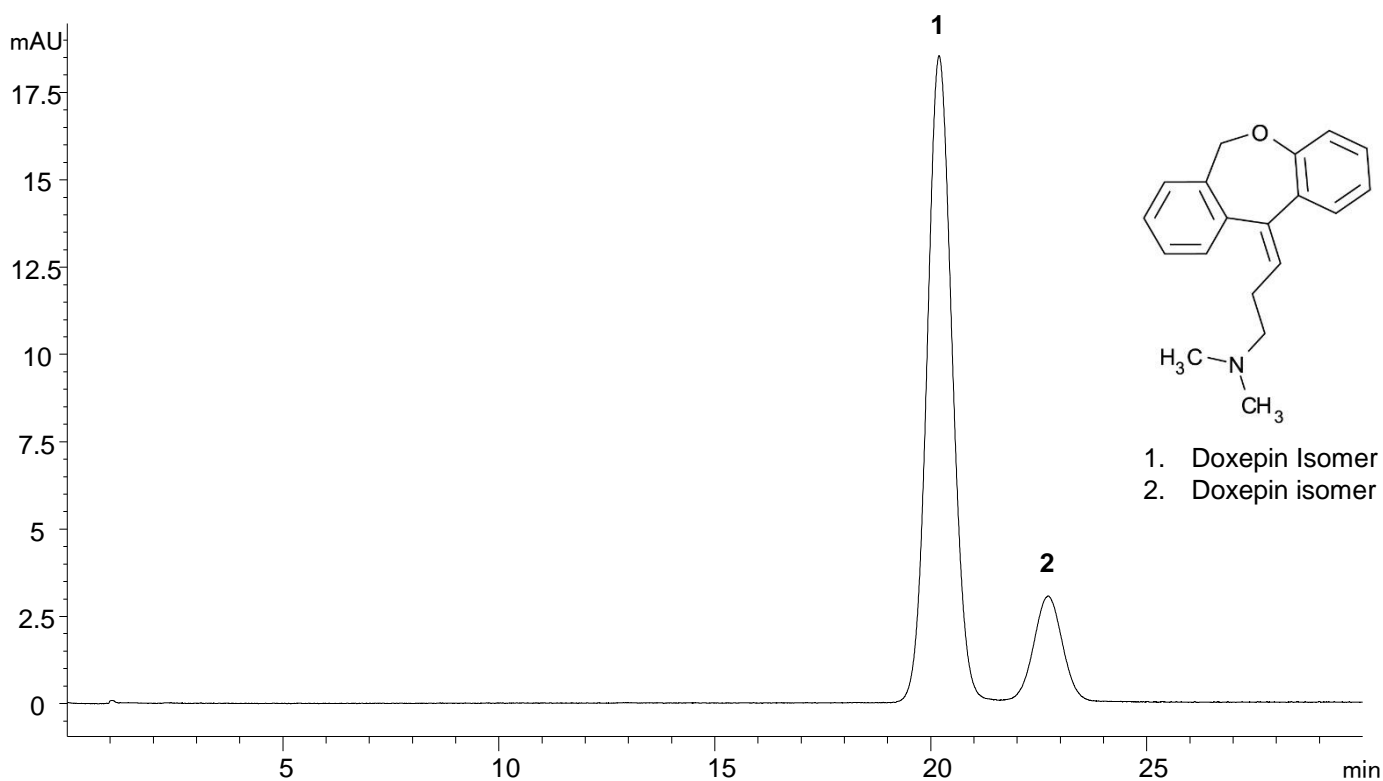

Translated method according to USP <621> guidance. Always check the latest and official method information from the relevant pharmacopoeia prior to analysis.

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.