

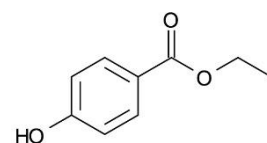


Conditions


Column:	ACE UltraCore 2.5 SuperC18
Dimensions:	100 x 4.6 mm
Part Number:	CORE-25A-1046U
Mobile Phase:	H ₂ O/MeCN (45:55 v/v)
Flow Rate:	1.39 mL/min
Injection:	10.1 µL
Temperature:	Ambient (22 °C)
Detection:	UV, 205 nm

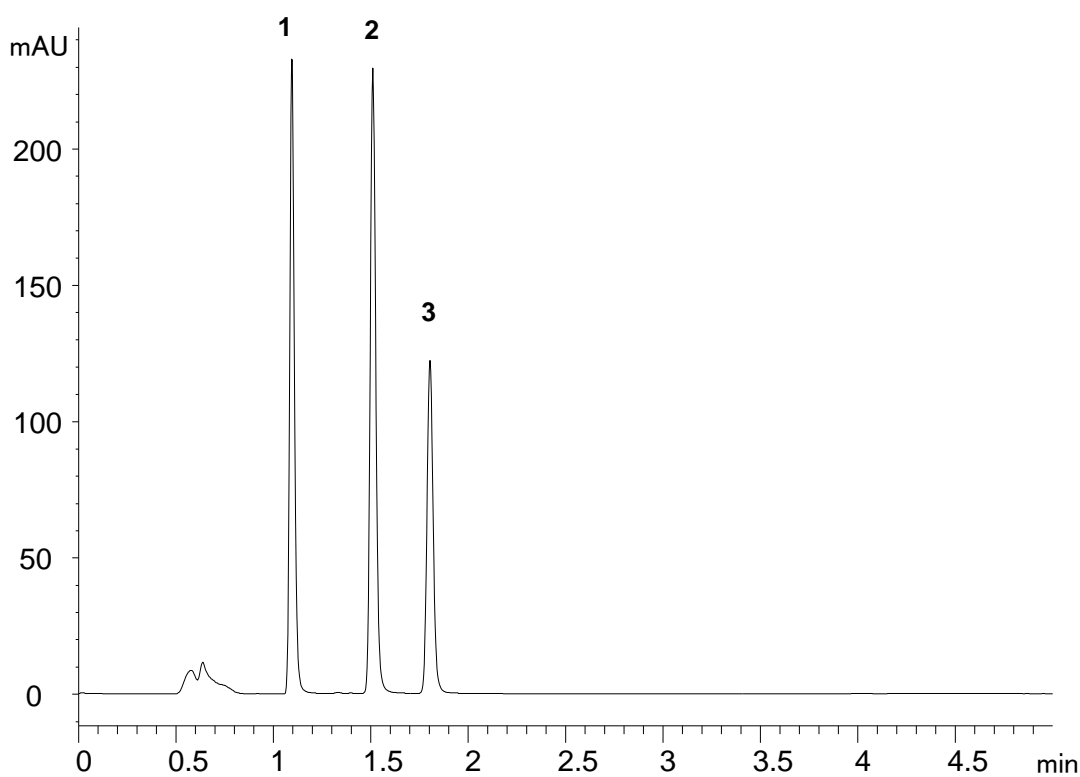

1. Ethylparaben

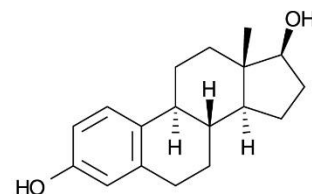

2. Estradiol

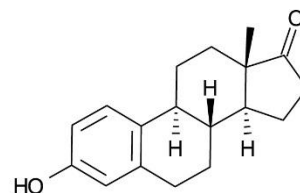

3. Estrone

Translated method according to USP <621> guidance. Always check the latest and official method information from the relevant pharmacopoeia prior to analysis.

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.