


Conditions

Column:	ACE 5 C18
Dimensions:	150 x 4.6 mm
Part Number:	ACE-121-1546
Mobile Phase:	8 mM monobasic sodium phosphate pH 2.4/MeCN (1:1 v/v)
Flow Rate:	1.32 mL/min
Injection:	16 µL
Temperature:	Ambient (22 °C)
Detection:	UV, 228 nm

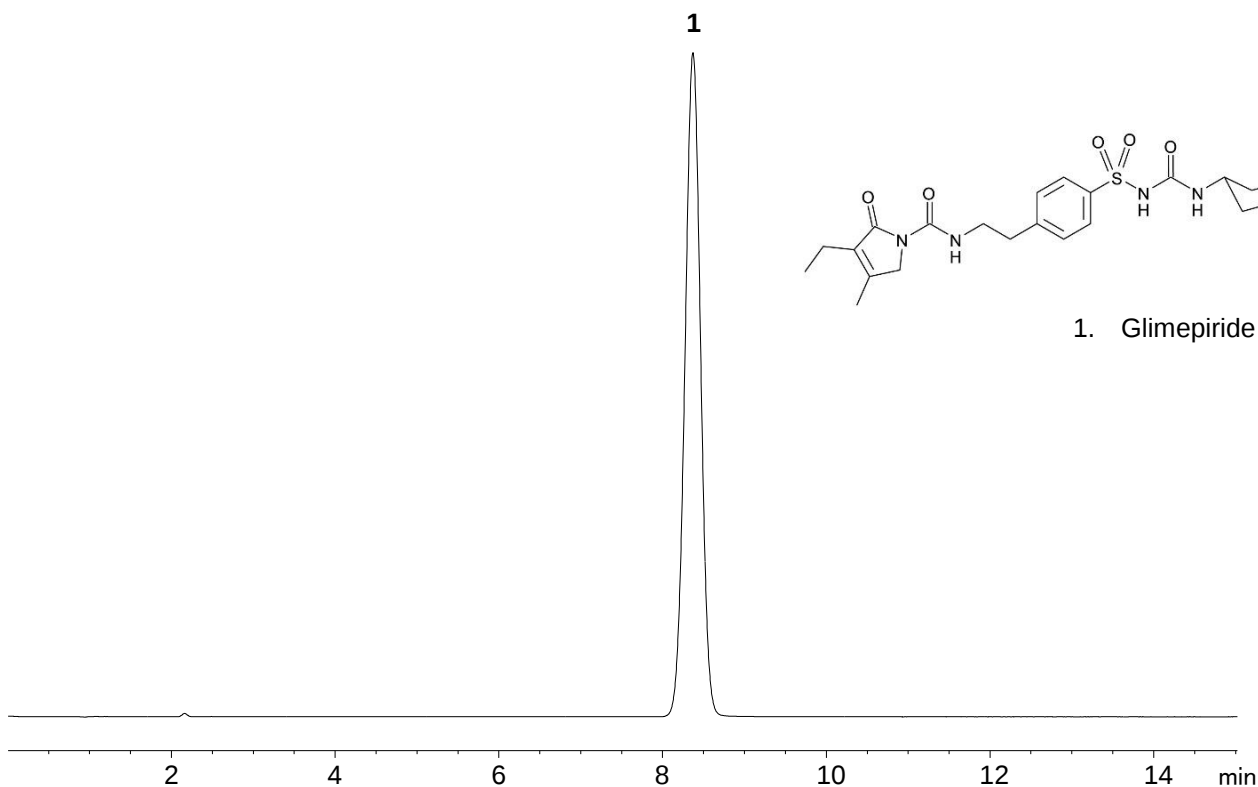

Translated method according to USP <621> guidance. Always check the latest and official method information from the relevant pharmacopoeia prior to analysis.

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com