



Vitamins – Fat Soluble

ACE®
Ultra-inert
UHPLC & HPLC Columns


Conditions

Column:	ACE 5 C18
Dimensions:	250 x 4.6 mm
Part Number:	ACE-121-2546
Mobile Phase:	MeOH
Flow Rate:	1.5 mL/min
Temperature:	30 °C
Detection:	UV, 280 nm


Application #AN6950

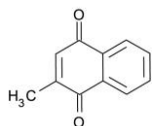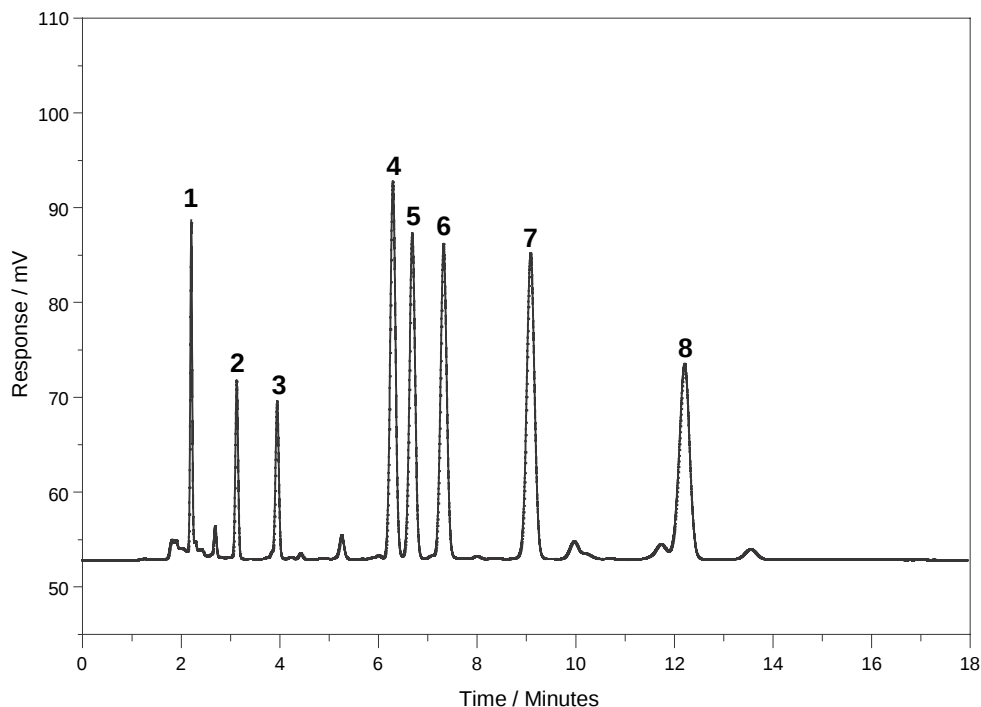


1. Menadione (Vitamin K3)

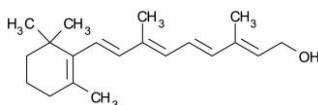


2. Retinol (Vitamin A)

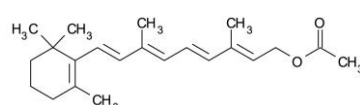


3. Vitamin A acetate

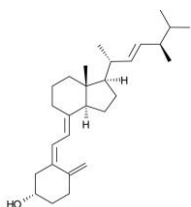


4. Ergocalciferol (Vitamin D2)

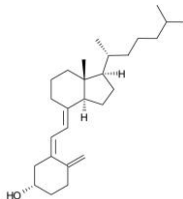


5. Cholecalciferol (Vitamin D3)

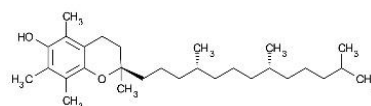

6. α -Tocopherol (Vitamin E)

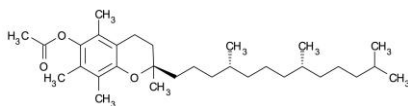

7. Vitamin E acetate

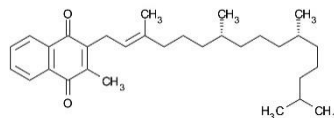

8. Vitamin K1

© 2019 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com



ACE
UHPLC & HPLC
Columns