




Vitamins – Fat Soluble

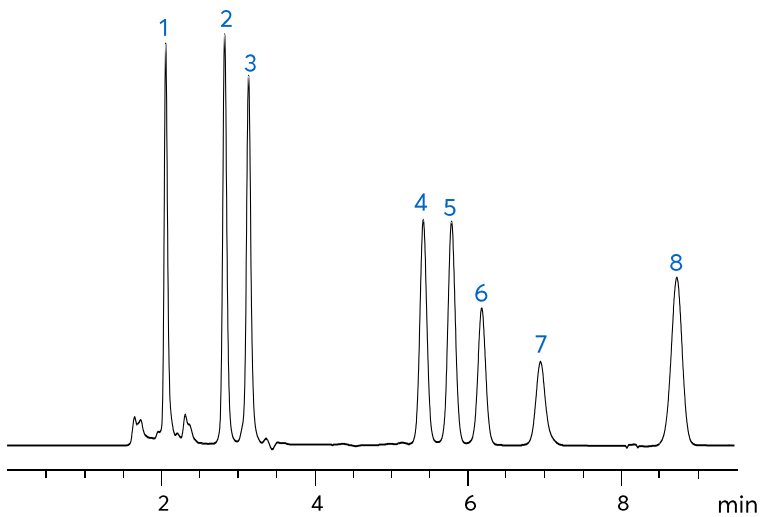


1. Menadione (Vitamin K3)
2. Retinol (Vitamin A)
3. Vitamin A acetate
4. Ergocalciferol (Vitamin D2)
5. Cholecalciferol (Vitamin D3)
6. Vitamin E acetate
7. α -Tocopherol (Vitamin E)
8. Vitamin K1

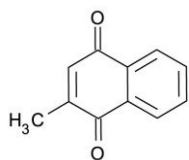


Menadione (Vitamin K3)

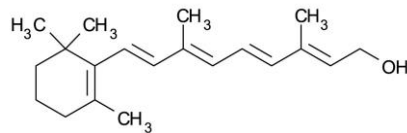


Retinol (Vitamin A)

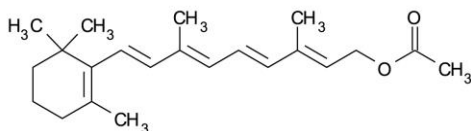


Vitamin A acetate

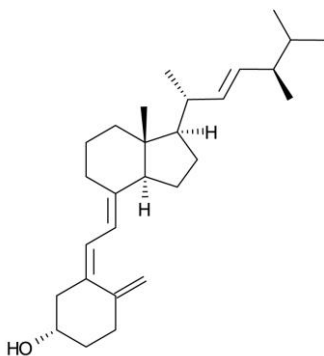


Ergocalciferol (Vitamin D2)

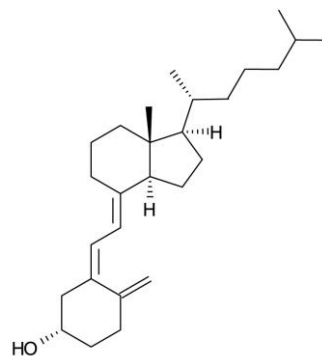


Cholecalciferol (Vitamin D3)

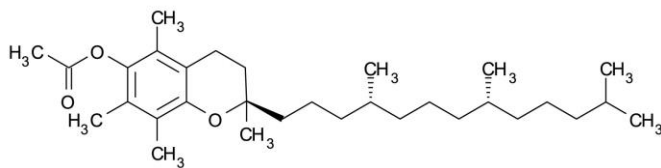

Vitamin E acetate

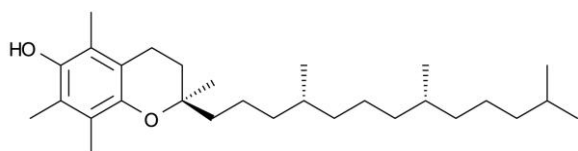

α-Tocopherol (Vitamin E)

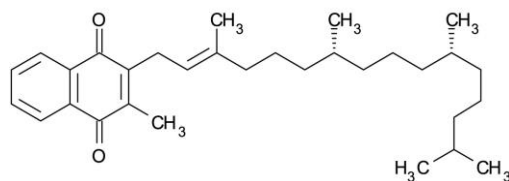

Vitamin K1

Method Details

CONDITIONS

Column:	Avantor® ACE® Excel® 3 C18-Amide
Particle Size:	3 µm
Dimensions:	150 x 4.6 mm
Mobile Phase:	MeOH/MeCN (90:10 v/v)
Flow Rate:	1 mL/min
Temperature:	20 °C
Detection:	UV, 280 nm

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® Excel® 3 C18-Amide	HPLC Column	150 x 4.6 mm	EXL-1112-1546U
Methanol	VWR CHROMANORM® gradient grade	2.5 L	20864.320
Acetonitrile	VWR CHROMANORM® super gradient grade	2.5 L	83639.320