


Vitamins – Water Soluble (Gradient 1)

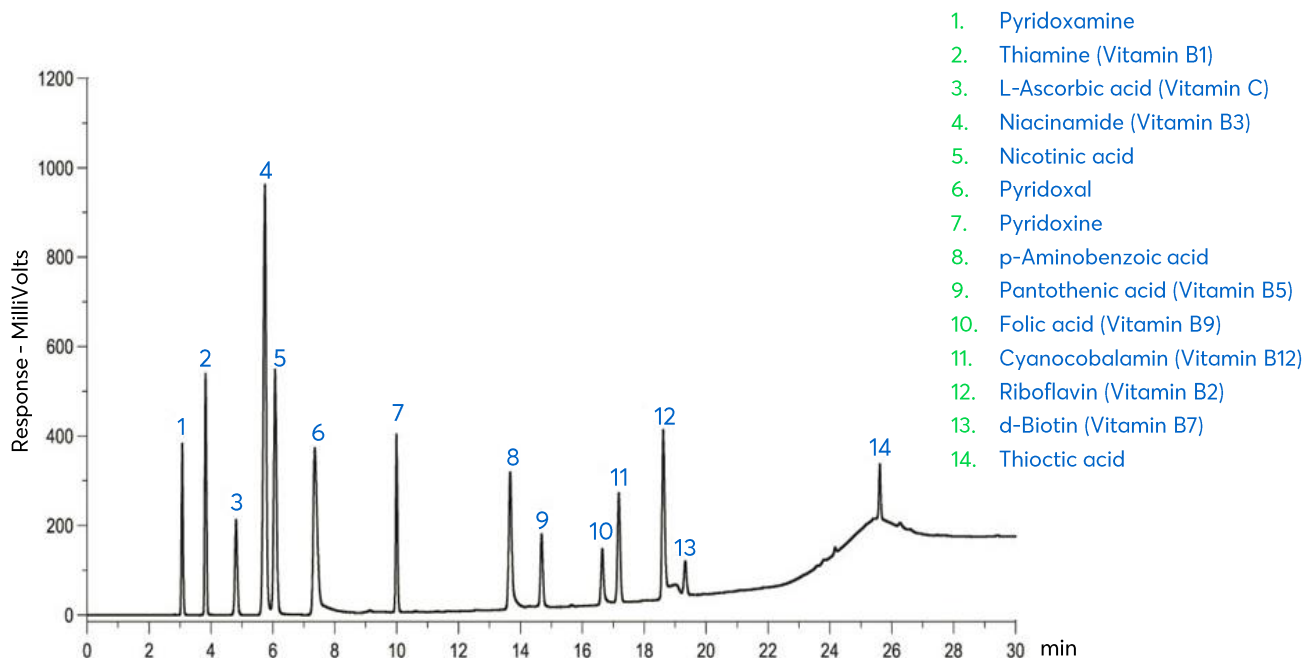

Method Details

CONDITIONS

Column: Avantor® ACE® 5 C8

Particle Size: 5 µm

Dimensions: 250 x 4.6 mm

Mobile Phase: A: 50 mM KH₂PO₄ pH 2.5 in H₂O
B: MeOH

Gradient:

Time (mins)	% B
0.0	0
3.0	0
16.5	45
19.5	80

Flow Rate: 1 mL/min

Temperature: Ambient

Detection: UV, 205 nm

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® 5 C8	HPLC Column	250 x 4.6 mm	ACE-122-2546
Methanol	VWR CHROMANORM® gradient grade	2.5 L	20864.320
Water	VWR CHROMANORM® super gradient grade	2.5 L	83650.320
Potassium dihydrogen phosphate	VWR CHROMANORM® for HPLC	500 g	153184U