


## Conditions

Column: ACE 3 C18  
Dimensions: 100 x 4.6 mm  
Part Number: ACE-111-1046  
Mobile Phase: A: 3 mM tetrabutylammonium bromide and 5 mM KH<sub>2</sub>PO<sub>4</sub> in H<sub>2</sub>O  
B: 5 mM tetrabutylammonium bromide in MeOH

Time (mins)	%B
0	45
20	70
30	45
40	45

Flow Rate: 0.8 mL/min  
Injection: 10 µL  
Temperature: Ambient  
Detection: UV-Vis, 420 nm, 520 nm and 600 nm

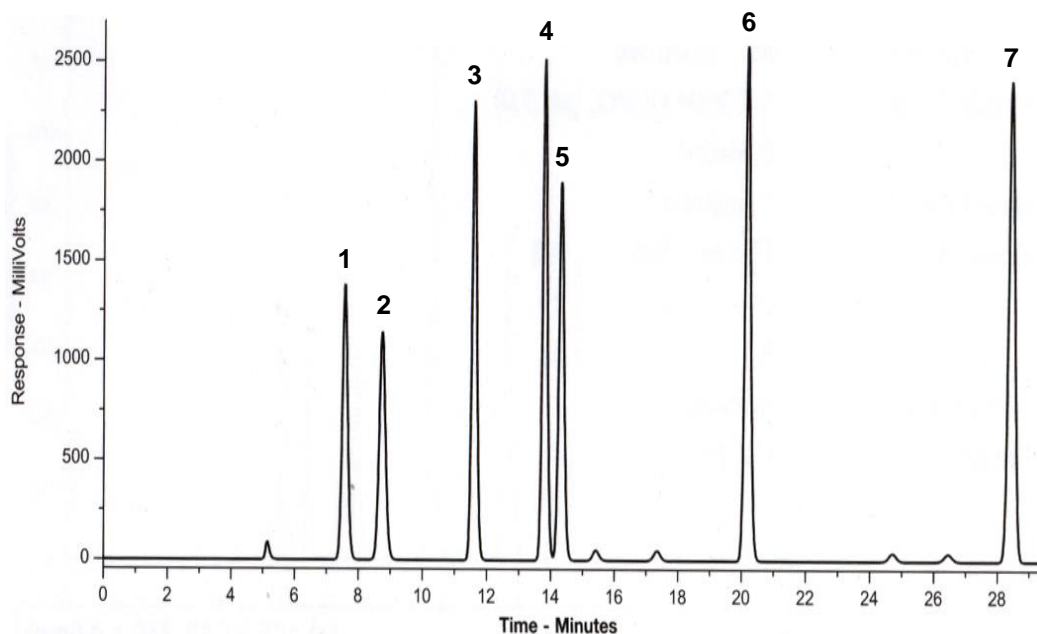

1. Amaranth
2. Sunset Yellow
3. Allura Red
4. Red 2G
5. Ponceau 4R
6. Carmoisine
7. Erythrosine