



mGluR5 PET Tracer by Radio HPLC Analysis

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN2700

Conditions

Column: ACE 3 C18
Dimensions: 50 x 4.6 mm
Part Number: ACE-111-0546
Mobile Phase: 0.1% TFA in H₂O/MeCN (55:45 v/v)
Flow Rate: 1 mL/min
Injection: 20 µL
Detection: UV, 254 nm
Radio detection

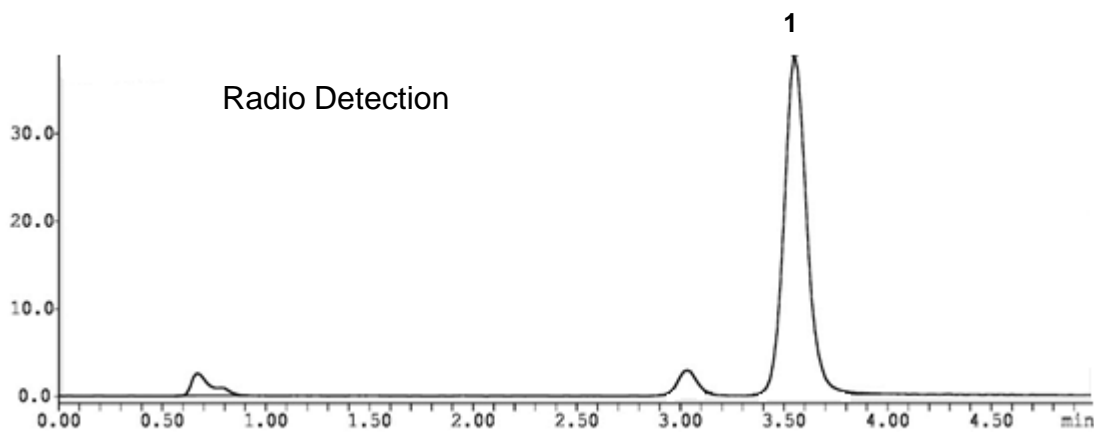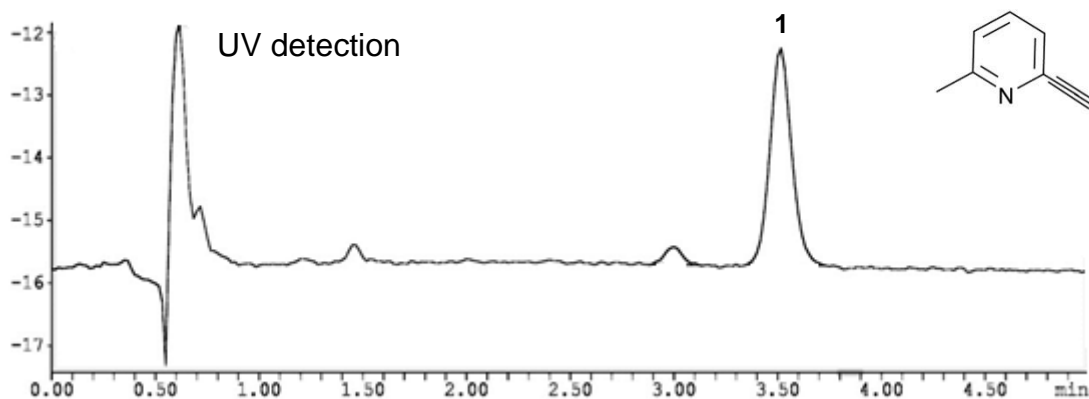

Reproduced with permission of Centre of Radiopharmaceutical Science, University Hospital Zurich, Switzerland

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.



For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com